

KANAMYCIN RAPID TEST STRIP
MANUAL

Kanamycin Rapid Test Strip (milk)

Catalogue Number. IP100098

Brief

The Kanamycin rapid test strip is used for testing Kanamycin residue in fresh milk qualitatively. The whole test only need 6min, east to operate with high sensitivity..

Principle

The Kanamycin rapid test strip is based on competitive inhibition immuno-chromatographic principle. In the flow process, Kanamycin in the sample combined with Kanamycin specific colloidal gold-labeled monoclonal antibody, Inhibit the combination between antibody and Kanamycin-BSA conjugate on Test line of NC membrane, lead to the color change of Test line. When the sample has no Kanamycin residue or concentration lower than detection limit, T line is darker then C line; when the concentration is equal to or higher than detection limit, T line is lighter than C line or T line has no color or T line has same color with C line. No matter whether there is Kanamycin residue in sample, C line will appear, it means the test is valid.

Technical specifications

8 strips/bottle, 12 bottles/kit

Detection limit: 3-5ppb

Components


1	Kanamycin test strip	8strip/bottle
2	Red powder micro-well	8 micro-well/bottle
3	Manual	1 piece
4	Desiccant	2 pieces/bag

Operation procedures

- 1) Read the instruction carefully before use. Return test dipsticks and sample into room temperature.
- 2) Take bottles needed from the kit package, take out required microwells, making proper marks. Seal the cap of the bottles, avoid moisture.
- 3) Take 200ul of the test samples into the microwells, then repeatedly absorb for 5 times, mix the sample with the reagent in the microwells completely until no solid judged by eyes (this is a very important step).
- 4) Incubate at 40±2°C for 3 min. Then insert the strip into micro-well, put the end with MAX down into the liquid.
- 5) Insert the strip into the liquid for 3min, then read the result.
- 7) Test Result Interpretation

Result Interpretation

Color depth comparison of T line and C line	Result judgment	Result analysis
T line > C line	Negative	Kanamycin residue is lower than detection limit in sample
T line = C line	Weak positive	Kanamycin residue is equal to detection limit in sample
T line < C line or T line is invisible in sample	Positive	Kanamycin residue is higher than detection limit in sample
C line is invisible	Invalid	The test dipstick is invalid, test again

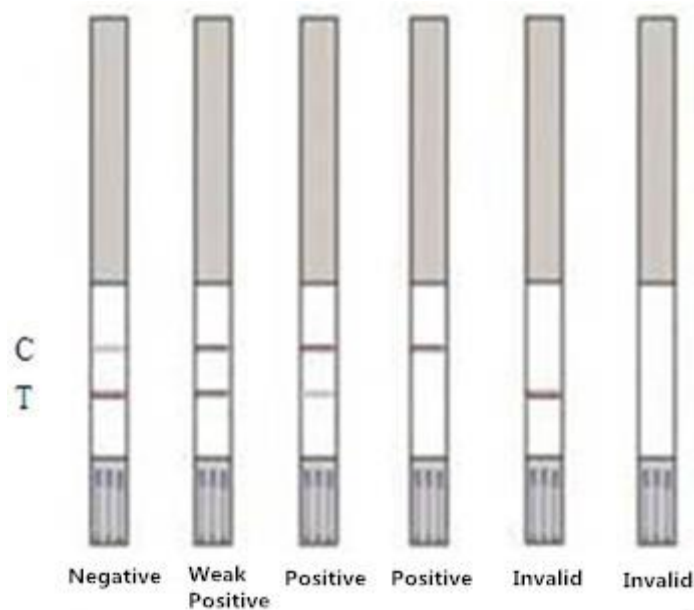

Precautions

- 1) The test strip is used for once at room temperature, do not use test strip out of date.
- 2) Cover the bottle immediately after taking out the reagent. If can't use 8 micro-wells for one time, cover the micro-wells immediately and put back to the bottle, store in seal.
- 3) Do not touch the white membrane surface in the middle of test strip, avoid sunlight and fan blowing directly.
- 4) Milk samples should be fully liquid without any agglomeration, sour and precipitation
- 5) Use the test strips testing again for positive results.
- 6) Please contact the supplier for any questions.
- 7) The result is for reference. To confirm, please refer national standard method.

Specificity

This product has no cross-reaction with Quinolones, Sulfonamides and Tetracyclines etc.

Storage: Store at 4-8 °C in dark, sealed, dry place, no frozen.

Expiry date: 12 months; date of production is on box..