


Conduritol B epoxide

Cat. No.:	HY-100944		
CAS No.:	6090-95-5		
Molecular Formula:	C ₆ H ₁₀ O ₅		
Molecular Weight:	162.14		
Target:	Glucosidase		
Pathway:	Metabolic Enzyme/Protease		
Storage:	Powder	-20°C	3 years
		4°C	2 years
	In solvent	-80°C	6 months
		-20°C	1 month

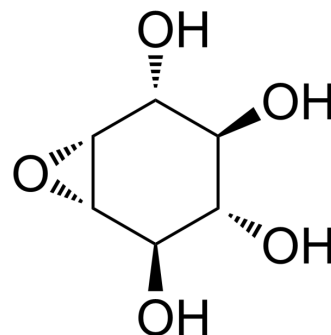

SOLVENT & SOLUBILITY

In Vitro

H₂O : 100 mg/mL (616.75 mM; Need ultrasonic)
 DMSO : 50 mg/mL (308.38 mM; Need ultrasonic)

	Solvent Concentration	Mass		
		1 mg	5 mg	10 mg
Preparing Stock Solutions	1 mM	6.1675 mL	30.8375 mL	61.6751 mL
	5 mM	1.2335 mL	6.1675 mL	12.3350 mL
	10 mM	0.6168 mL	3.0838 mL	6.1675 mL

Please refer to the solubility information to select the appropriate solvent.

In Vivo

- Add each solvent one by one: PBS
Solubility: 100 mg/mL (616.75 mM); Clear solution; Need ultrasonic
- Add each solvent one by one: 10% DMSO >> 40% PEG300 >> 5% Tween-80 >> 45% saline
Solubility: ≥ 2.5 mg/mL (15.42 mM); Clear solution
- Add each solvent one by one: 10% DMSO >> 90% (20% SBE-β-CD in saline)
Solubility: ≥ 2.5 mg/mL (15.42 mM); Clear solution
- Add each solvent one by one: 10% DMSO >> 90% corn oil
Solubility: ≥ 2.5 mg/mL (15.42 mM); Clear solution

BIOLOGICAL ACTIVITY

Description

Conduritol B epoxide is an irreversible covalently bound acid β-glucosidase (GCCase) inhibitor.

In Vitro

Treating N2a and Conduritol B epoxide-N2a cells with G6 gave GC reductions to 7% and 26% of untreated levels, respectively. G6 treated Conduritol B epoxide-N2a also has significantly decreased GS. Co-treatment of G6 and dantrolene of

Conduritol B epoxide-N2a cells lead to a similar degree of GC and GS reduction as G6 alone, indicating a specific effect of G6 on inhibition of substrate accumulation, Conduritol B epoxide-N2a cells have higher calcium levels than in N2a cells without caffeine at baseline. Conduritol B epoxide-N2a cells show a significant increase in cytosolic calcium levels compared to N2a cells. Conduritol B epoxide-N2a cells show a significant reduction in OCR as evidenced by approximately 50% in all the parameters, including rate of ATP production, basal respiration, and maximal respiration, compared to N2a cells, indicating reduced mitochondrial function in this nGD cell model. In dantrolene-treated Conduritol B epoxide-N2a cells, levels of Ryr3 are increased to 76% of WT level^[1].

MCE has not independently confirmed the accuracy of these methods. They are for reference only.

In Vivo

Conduritol B epoxide treatment of 4L, 9H, 9V and WT mice (100 mg/kg/day, i.p.) from postnatal day 5 to 11 does not induce α -synuclein aggregates. Long-term daily Conduritol B epoxide treatment of 4L mice (100 mg/kg/day, from postnatal day 15 for 24 or 36 doses) leads to hind limb paralysis and small amounts of α -synuclein accumulation in the olfactory bulb, brainstem, and PVP near D3V (dorsal 3rd ventricle)^[2].

MCE has not independently confirmed the accuracy of these methods. They are for reference only.

PROTOCOL

Kinase Assay ^[1]

Cells are homogenized in 1% sodium taurocholate/1% Triton X-100. GCCase activity is determined fluorometrically using 4MU-Glucose as the substrate in 0.25% sodium taurocholate, 0.25% Triton X-100 and 0.1M citric-phosphate buffer (pH 5.6). Brain tissues are homogenized in 1X PBS and incubated in 5 μ M brain phosphatidylserine and 0.1M citric-phosphate buffer (pH 5.6) for GCCase activity assay using 4MU-Glucose as substrate. Protein concentrations are determined by BCA assay using BSA as standard.

MCE has not independently confirmed the accuracy of these methods. They are for reference only.

Animal Administration ^[2]

Gba1 point mutated 4L, 9H and 9V mice or WT mice are intraperitoneally injected with 100 mg/kg/day of Conduritol B epoxide. In short-term experiments, daily injections are initiated at postnatal day 5 and continued for 6 daily doses. The mice are sacrificed on day 12 or 2 months after last injection. In long-term experiments, 4L mice are injected daily beginning at postnatal day 15 for 24 or 36 daily doses and sacrificed the day after last injection. Mice are perfused with PBS and organs are harvested for enzyme activity, lipid and histological analyses.

MCE has not independently confirmed the accuracy of these methods. They are for reference only.

CUSTOMER VALIDATION

- Neuron. 2020 Mar 4;105(5):822-836.e7.
- Mol Biol Evol. 2022 Feb 22;msac042.
- Inflammation. 2021 Mar 4.
- Toxins. 2022, 14(7), 458.
- RSC Adv. 2018 8:8469-8483.

See more customer validations on www.MedChemExpress.com

REFERENCES

[1]. Liou B, et al. Modulating ryanodine receptors with dantrolene attenuates neuronopathic phenotype in Gaucher disease mice. Hum Mol Genet. 2016 Sep 20. pii: ddw322. [Epub ahead of print]

[2]. YH Xu, et al. Accumulation and distribution of α -synuclein and ubiquitin in the CNS of Gaucher disease mouse models. Mol Genet Metab. Author manuscript; available

Caution: Product has not been fully validated for medical applications. For research use only.

Tel: 609-228-6898

Fax: 609-228-5909

E-mail: tech@MedChemExpress.com

Address: 1 Deer Park Dr, Suite Q, Monmouth Junction, NJ 08852, USA