


AZD2858

Cat. No.:	HY-15761		
CAS No.:	486424-20-8		
Molecular Formula:	C ₂₁ H ₂₃ N ₇ O ₃ S		
Molecular Weight:	453.52		
Target:	GSK-3		
Pathway:	PI3K/Akt/mTOR; Stem Cell/Wnt		
Storage:	Powder	-20°C	3 years
		4°C	2 years
	In solvent	-80°C	1 year
		-20°C	6 months

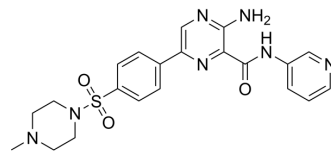

SOLVENT & SOLUBILITY

In Vitro	DMSO : 12.5 mg/mL (27.56 mM; Need ultrasonic)				
		Solvent Concentration	Mass 1 mg	5 mg	10 mg
	Preparing Stock Solutions	1 mM	2.2050 mL	11.0249 mL	22.0497 mL
		5 mM	0.4410 mL	2.2050 mL	4.4099 mL
10 mM		0.2205 mL	1.1025 mL	2.2050 mL	
Please refer to the solubility information to select the appropriate solvent.					
In Vivo	1. Add each solvent one by one: 10% DMSO >> 40% PEG300 >> 5% Tween-80 >> 45% saline Solubility: 1.25 mg/mL (2.76 mM); Suspended solution; Need ultrasonic 2. Add each solvent one by one: 10% DMSO >> 90% (20% SBE-β-CD in saline) Solubility: ≥ 1.25 mg/mL (2.76 mM); Clear solution				

BIOLOGICAL ACTIVITY

Description	AZD2858 is a potent, orally active GSK-3 inhibitor, with IC ₅₀ s of 0.9 and 5 nM for GSK-3α and GSK-3β, respectively, used in the research of fracture healing.			
IC₅₀ & Target	GSK-3α 0.9 nM (IC ₅₀)	GSK-3β 5 nM (IC ₅₀)	CDK5/p25 356 nM (IC ₅₀)	Haspin 366 nM (IC ₅₀)
	CDK5/p35 387 nM (IC ₅₀)	DYRK2 491 nM (IC ₅₀)	CDK2/cyclin A 810 nM (IC ₅₀)	CDK1/cyclin B 1246 nM (IC ₅₀)
	PIM3	TLK2	PKD2	CDK2/cyclin E

	1269 nM (IC ₅₀)	1381 nM (IC ₅₀)	2462 nM (IC ₅₀)	3310 nM (IC ₅₀)
	Aurora-A 4966 nM (IC ₅₀)			
In Vitro	<p>AZD2858 (1 μM) increases β-catenin levels after a short period of time in human osteoblast cells. AZD2858 inhibits GSK-3β dependent phosphorylation with an IC₅₀ of 68 nM. AZD2858 (10 nM) has no effect on β-catenin levels^[1]. AZD2858 increases TAZ expression and osterix expression both by 1.4-fold, with EC₅₀ of 440 nM and 1.2 μM, respectively, in hADSC. AZD2858 also induces a marked increase in osteogenic mineralisation in hADSC^[3]. AZD2858 (AR28) demonstrates from 70- to greater than 6000-fold selectivity over a panel of other kinases and an IC₅₀ of 5 nM. AR28 inhibits GSK-3 in murine cells and indicates activation of the canonical Wnt/β-catenin signaling cascade. AR28 (50, 10, and 1 nM) enhances the clonogenic ability of mesenchymal progenitors with osteogenic and adipogenic potential. AR28 (50 μM) also enhances the differentiation ability of mesenchymal progenitors to the osteogenic but not adipogenic lineage in vitro^[4].</p> <p>MCE has not independently confirmed the accuracy of these methods. They are for reference only.</p>			
In Vivo	<p>AZD2858 (20 mg/kg) causes a dose-dependent increase in trabecular bone mass compared to control after a two-week treatment with a maximum effect^[1]. AZD2858 exhibits a substantial effect on fracture healing. AZD2858 (20 mg/kg) causes an increase in cortical BMC of 9%, cortical area of 10%, and cortical thickness of 11% at 3 weeks in the non-operated right femur of rats^[2]. AZD2858 (30 μmol/kg/day) alters the biomarkers of bone turnover with statistically significant increases in P1NP and decreases in TRAcP-5b seen from 3 days of treatment and onwards. AZD2858 demonstrates significant changes in serum bone turnover markers (P1NP and TRAcP-5b) and femur bone formation after only 7 days of daily dosing^[3]. AZD2858 (AR28, 30 mg/kg, s.c.) stimulates an increase in an initial wave of mesenchymal progenitors with osteogenic and adipogenic potential and drives their differentiation to the osteogenic lineage in BALB/c mice. AR28 (30 mg/kg, s.c.) enhances the proliferation of committed hematopoietic progenitors and their differentiation to the osteoclast lineage but does not prevent an overall increase in bone mass^[4].</p> <p>MCE has not independently confirmed the accuracy of these methods. They are for reference only.</p>			

PROTOCOL

Kinase Assay ^[3]

The potency of compounds at GSK-3β and cyclin-dependent protein kinase 2 (CDK2, kinase with closest homology to GSK-3 β) is assessed using Z-LYTE™ Kinase assay kit in the presence of 7 and 80 μM ATP respectively. A ratiometric method is used to calculate the ratio of donor emission (445 nm) to acceptor emission (520 nm) after excitation of the donor fluorophore at 400 nm to quantitate the reaction progress. Kinase selectivity with AR79, AZD2858 and AZ13282107 are determined using the KinaseProfiler Service or University of Dundee Kinase. Over 80 different kinases are assessed at a single concentration of 1 or 10 μM of AR79, AZD2858 and AZ13282107. Concentration-inhibition 10-point curves to compounds that show activity are constructed to determine pIC₅₀ estimations. Also, in some kinase assays these pIC₅₀ estimations are converted to binding affinity values (pK_i) using the Cheng-Prussoff equation to correct for the concentration of ATP used^[3].

MCE has not independently confirmed the accuracy of these methods. They are for reference only.

Animal Administration ^[1]

Each rat is dosed orally with vehicle or AZD2858 using a plastic gavage tube. The dose volume is 10 mL/kg. The vehicle consists of deionized water adjusted to pH 3.5±0.1. Formulations are adjusted to pH 3.5±0.1. The doses are 0, 0.2, 2 or 20 mg/kg respectively and administered either twice daily (TD), once daily (OD), every other day (O/2D), or every fourth day (O/4D) for 14 days. In total, the protocol results in 13 groups with 8 animals in each group (104 animals). At 7 days after the start of the study, and again 2 days prior to the scheduled terminal necropsy, each animal is injected subcutaneously with a bicarbonate buffered calcein solution (8 mg/kg, 1 mL/kg)^[1].

MCE has not independently confirmed the accuracy of these methods. They are for reference only.

CUSTOMER VALIDATION

- Cell Res. 2022 Jun;32(6):513-529.
- Oncogene. 2023 Oct 16.
- J Cell Physiol. 2019 Jan 26.
- J Biol Chem. 2016 Dec 2;291(49):25529-25541.
- Oncol Rep. 2019 Jul;42(1):91-102.

See more customer validations on www.MedChemExpress.com

REFERENCES

- [1]. Marsell R, et al. GSK-3 inhibition by an orally active small molecule increases bone mass in rats. Bone. 2012 Mar;50(3):619-27.
- [2]. Sisask G, et al. Rats treated with AZD2858, a GSK3 inhibitor, heal fractures rapidly without endochondral bone formation. Bone. 2013 May;54(1):126-32.
- [3]. Gilmour PS, et al. Human stem cell osteoblastogenesis mediated by novel glycogen synthase kinase 3 inhibitors induces bone formation and a unique bone turnover biomarker profile in rats. Toxicol Appl Pharmacol. 2013 Oct 15;272(2):399-407.
- [4]. Gambardella A, et al. Glycogen synthase kinase-3 α/β inhibition promotes in vivo amplification of endogenous mesenchymal progenitors with osteogenic and adipogenic potential and their differentiation to the osteogenic lineage. J Bone Miner Res. 2011 Apr;26(4):811-21.
-

Caution: Product has not been fully validated for medical applications. For research use only.

Tel: 609-228-6898

Fax: 609-228-5909

E-mail: tech@MedChemExpress.com

Address: 1 Deer Park Dr, Suite Q, Monmouth Junction, NJ 08852, USA