


Voreloxin Hydrochloride

Cat. No.:	HY-16518
CAS No.:	175519-16-1
Molecular Formula:	C ₁₈ H ₂₀ ClN ₅ O ₄ S
Molecular Weight:	437.9
Target:	Topoisomerase; Apoptosis
Pathway:	Cell Cycle/DNA Damage; Apoptosis
Storage:	4°C, sealed storage, away from moisture * In solvent : -80°C, 6 months; -20°C, 1 month (sealed storage, away from moisture)

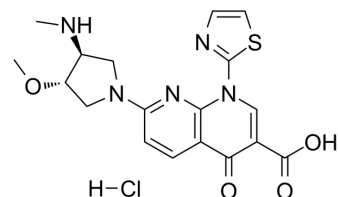

SOLVENT & SOLUBILITY

In Vitro	DMSO : 2 mg/mL (4.57 mM); ultrasonic and warming and heat to 60°C																					
	H ₂ O : 2 mg/mL (4.57 mM); ultrasonic and warming and heat to 60°C																					
Preparing Stock Solutions	<table border="1"> <thead> <tr> <th rowspan="2">Solvent</th> <th rowspan="2">Mass</th> <th>1 mg</th> <th>5 mg</th> <th>10 mg</th> </tr> <tr> <th>Concentration</th> <th></th> <th></th> <th></th> </tr> </thead> <tbody> <tr> <td>1 mM</td> <td>2.2836 mL</td> <td>11.4181 mL</td> <td>22.8363 mL</td> </tr> <tr> <td>5 mM</td> <td>---</td> <td>---</td> <td>---</td> </tr> <tr> <td>10 mM</td> <td>---</td> <td>---</td> <td>---</td> </tr> </tbody> </table>	Solvent	Mass	1 mg	5 mg	10 mg	Concentration				1 mM	2.2836 mL	11.4181 mL	22.8363 mL	5 mM	---	---	---	10 mM	---	---	---
	Solvent			Mass	1 mg	5 mg	10 mg															
		Concentration																				
	1 mM	2.2836 mL	11.4181 mL	22.8363 mL																		
5 mM	---	---	---																			
10 mM	---	---	---																			
Please refer to the solubility information to select the appropriate solvent.																						
In Vivo	1. Add each solvent one by one: PBS Solubility: 2 mg/mL (4.57 mM); Clear solution; Need ultrasonic and warming and heat to 60°C																					

BIOLOGICAL ACTIVITY

Description	Voreloxin Hydrochloride is a first-in-class topoisomerase II inhibitor that intercalates DNA and induces site-selective DNA DSB, G2 arrest, and apoptosis.
IC₅₀ & Target	Topoisomerase II
In Vitro	Voreloxin Hydrochloride is a first-in-class topoisomerase II poison and inhibitor that intercalates DNA and induces site-selective DNA DSB, G2 arrest, and apoptosis. Voreloxin (0.1-20 μM) inhibits topoisomerase II activity and induces site-selective DNA DSB in CCRF-CEM cells. Voreloxin (0.11, 0.33, 1, 3 μM) induces G2 arrest partially through topoisomerase II in A549 lung cancer cell line. Voreloxin cytotoxic activity requires DNA intercalation. However, Voreloxin (1-9 μM) does not generate significant levels of ROS ^[1] . Voreloxin has potent cytotoxic activity in AML cell lines MV4-11 and HL-60, with IC ₅₀ s of 95 ± 8 nM and 884 ± 114 nM, respectively. Voreloxin in combination with cytarabine shows additive or synergistic activity in acute leukemia cell lines ^[2] . Voreloxin is active on the primary acute myeloid leukemia (AML) with a mean LD ₅₀ of 2.3 μM. The LD ₅₀ for voreloxin in myeloid cell lines NB4 and HL-60 is 0.59 μM ± 0.25 μM. Voreloxin causes accumulation of cells in the S

and G2 phases of the cell cycle and acts on topoisomerase II^[3].

MCE has not independently confirmed the accuracy of these methods. They are for reference only.

In Vivo

Voreloxin (20 mg/kg, i.v.) alone results in 80% reduction in bone marrow cellularity of CD-1 mice by administration one dose every 4 days repeated twice (q4d ×2). voreloxin at 10 mg/kg in combination with cytarabine causes ablation of the marrow, dilation of sinusoids, and infiltration of adipocytes in mice. Voreloxin (20 mg/kg, i.v.) combined with cytarabine causes a reversible decrease in myeloid and lymphoid cells in bone marrow and peripheral blood CD-1 mice. voreloxin (10 mg/kg, q4d ×2) and cytarabine in combination causes reversible neutropenia with a more modest impact on platelets CD-1 mice^[2]. MCE has not independently confirmed the accuracy of these methods. They are for reference only.

PROTOCOL

Cell Assay ^[3]

In vitro toxicity assays are performed on primary AML mononuclear cells over a 48 h period using a MTS cell proliferation assay. Lethal doses (LD₅₀) are calculated. Cells are treated with voreloxin (31.25 nM to 4 μM) and Ara-C (62.5 nM to 8 μM) by serial dilution and incubated for 48 h in a final volume of 90 μL. Following the incubation, 20 μL of MTS reagent are added and the reaction is incubated for a further 4 h. The absorbance of the reaction after this time is read by spectrophotometry at 490 nm and the percentage of viable cells calculated relative to untreated control cells in the same assay. IC₅₀ values are calculated using Calcsyn software^[3].

MCE has not independently confirmed the accuracy of these methods. They are for reference only.

Animal Administration ^[2]

Animals are weighed, randomized by body weight, and assigned to the study groups before initiation of treatment. Voreloxin is administered intravenously (IV) at 10 or 20 mg/kg once on day zero and once on day four (q4d ×2). Cytarabine is administered subcutaneously (SC) at 20 or 60 mg/kg every 8 h on day zero and day four (tid q4d ×2). Tissues and blood are sampled on days 6, 8, and 12 from at least three and not greater than ten animals per treatment group. Femurs are placed in Streck Tissue Fixative solution, or in 10% formalin solution for 24-48 h followed by a 70% dehydrant (ethanol, isopropanol, methanol). Femurs are decalcified, paraffin embedded, and sectioned at Biopathology Labs. The four micron sections are stained with hematoxylin-eosin (H&E). H&E stained femurs are examined and percent cellularity of the bone marrow is determined. Digital photographs of representative femur sections are taken on a Leica DM2000 microscope using Image-Pro Plus v6.1 software^[2].

MCE has not independently confirmed the accuracy of these methods. They are for reference only.

CUSTOMER VALIDATION

- Drug Des Devel Ther. 2018 Apr 30;12:1009-1017.
- Oncotarget. 2018 Mar 30;9(24):16917-16931.

See more customer validations on www.MedChemExpress.com

REFERENCES

[1]. Hotinski AK, et al. Vosaroxin is a novel topoisomerase-II inhibitor with efficacy in relapsed and refractory acute myeloid leukaemia. Expert Opin Pharmacother. 2015 Jun;16(9):1395-402.

[2]. Scatena CD, et al. Voreloxin, a first-in-class anticancer quinolone derivative, acts synergistically with cytarabine in vitro and induces bone marrow aplasia in vivo. Cancer Chemother Pharmacol. 2010 Oct;66(5):881-8.

[3]. Walsby EJ, et al. The topoisomerase II inhibitor voreloxin causes cell cycle arrest and apoptosis in myeloid leukemia cells and acts in synergy with cytarabine. Haematologica. 2011 Mar;96(3):393-9.

Caution: Product has not been fully validated for medical applications. For research use only.

Tel: 609-228-6898

Fax: 609-228-5909

E-mail: tech@MedChemExpress.com

Address: 1 Deer Park Dr, Suite Q, Monmouth Junction, NJ 08852, USA